



MEDICONS Cod. 79300N / 79300S

Disposable disaggregator with 35µm cut microblades stainless steel disk. HD/LD PE casing. Provided sterile and non-sterile.

The stainless steel disk is made up of a fixed screen with about 100 hexagonal bore-holes, around each of which there are 6 microblades.

- 1 ml liquid solution
- All kind of fresh/frozen solid tissues can be processed

USE

The non sterile Medicons code 79300N can be used for general laboratory use and for research to obtain live single cells of size below 50µm.

The sterile Medicons code 79300S can be used for general laboratory use, research and diagnostics to obtain live single cells of size below 50µm.

Its stainless steel micro-blades positioned at hexagonal pyramid around each hole of the grinder disk provide a cut size of about 35µm; able to disaggregate all kind of tissues (fresh/frozen) including bone and cartilage samples.

Sample size: small parts of about 2 mm³


Dilution: 1 to 1,2 ml

Common tissues: Skin, Lungs, Kidneys, Bones, Cartilage

The Medicons 79300N can be run with the below devices:



RUO

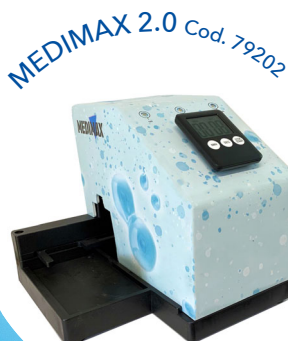

IVD



MD

www.ctsv.biz

